
truelife

SonicBrush Clean70 UV

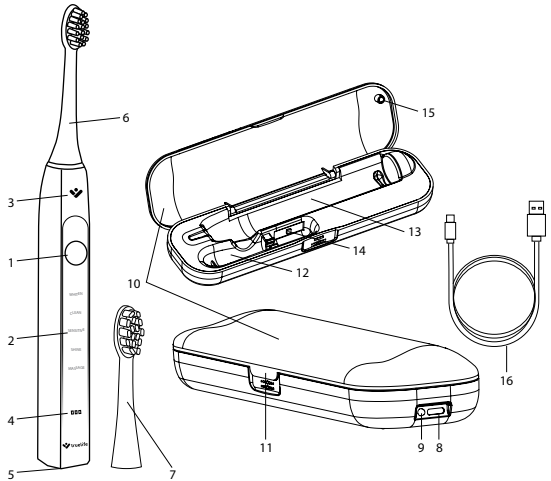
EN – INSTRUCTIONS

- PACKAGE CONTENTS**
- TrueLife SonicBrush Clean70 UV toothbrush
 - Rechargeable travel case with UVC
 - 1x Standard head, 1x Sensitive head
 - USB-C cable
 - Instructions
 - Safety instructions

- PRODUCT DESCRIPTION**
- ON / OFF button and mode change
 - Modes indicator
 - Pressure sensor
 - Battery indicator
 - USB-C input for power cable
 - Sensitive head
 - Head standard
 - USB-C input for power cable
 - Battery/sterilization indicator
 - Travel case lid
 - Case lock
 - Head hole
 - Hole for brush
 - UVC light
 - Magnet
 - USB-C cable

TEETH CLEANING

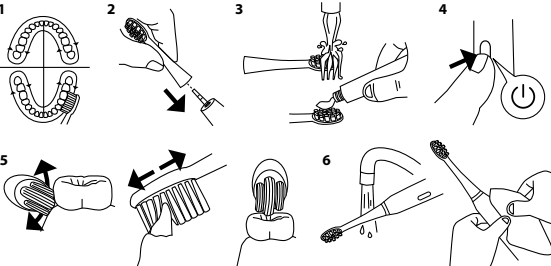- For maximum efficiency, divide the cleaning into 4 sections and devote 30 seconds to each of them. After 30 seconds, the brush will take a short pause and alert you to move to the next section. After 2 minutes, the brush turns itself off. Cleanse at least twice a day, in the morning and in the evening.
- Fit the brush head firmly onto the metal shaft with the bristles facing the front of the handle.
- Moisten the head, apply toothpaste and rub it on the teeth.
- Switch on the brush with the ON / OFF button. To select another mode, press the button again and select one of the five modes. If you would like to end cleaning earlier, hold the ON/OFF button.
- In each section, start by cleaning the outside of the teeth, then the inside of the teeth, and finally clean the surfaces. Hold the brush at an angle so that the bristles gently touch the teeth and the edges of the gums at the same time. When cleaning, do not press the brush and do not make circular movements like with a manual brush, but move slowly from one tooth to another.
- Rinse the brush head and body under running water and then wipe them with a cloth or towel.

For more information on how to brush your teeth properly:

truelife.eu/01921/video

CLEANING MODES

- WHITEN** – For removing stubborn dental plaque and whitening teeth
CLEAN – For routine cleaning and plaque removal
SENSITIVE – For cleaning sensitive teeth and gums
SHINE – To remove stains from coffee and tea and achieve a brilliant shine of teeth
MASSAGE – To stimulate the gums and improve their health

REPLACEMENT OF HEADS

We recommend changing the heads every 3 months. After this period of use, wear occurs where the frayed and bent fibers have a reduced ability to remove dental plaque and bacteria accumulate on them.

PRESSURE SENSOR

The brush has a pressure sensor on the front and back of the brush.
Blue light: suitable cleaning pressure
Green light: maximum usable pressure in a healthy range
Red light: excessive pressure (the brush warns with intermittent vibrations)

Note: the intermittent vibration when too much pressure can be turned off by holding the ON / OFF button for 5 seconds when the brush is turned off.

CHARGING THE TOOTHBRUSH

- When the brush is almost discharged, the charging indicator lights up red. When the brush is fully discharged, the charging indicator flashes red and cannot be used.
- Connect the connector of the USB cable to the socket of the brush, the other end to the adapter and the adapter to the outlet. You can also charge the brush from a computer or power bank. Charge the sonic brush using the included USB-C cable. Do not use fast chargers for charging.
- During charging, the indicator is loaded in white. Once the brush is charged, the charging indicator lights up white. Charging time is 4-5 hours.

CASE WITH UVC LIGHT

Insert the head into the housing with the bristles facing the UVC light. After closing the lid, the sterilization process begins, which lasts exactly 5 minutes. During sterilization, the indicator lights up green. When sterilization is complete, the indicator will turn off. Do not open the case lid during sterilization.

CHARGING THE CASE

- When the case is almost empty, the charging indicator will light red. When the case is completely discharged, the charge indicator flashes red and needs to be charged.
- Connect the connector of the USB cable to the socket of the brush, the other end to the adapter and the adapter to the socket. You can also charge the brush from a computer or power bank. Charge the sonic brush using the included USB-C cable. Do not use fast chargers for charging.

- During charging, the indicator lights red. When the case is charged, the charging indicator lights green. Charging time is 2-3 hours.

Note: the charging case is not used to charge the brush, but only to sterilize the heads.

TOOTHBRUSH SPECIFICATIONS

Battery life	70 days (2x daily 2 minutes)
Charging	USB-C cable
Number of moves	90,000 movements / min
Modes	Whiten, Clean, Sensitive, Shine, Massage
Pressure sensor	Yes
Timer	2 min (30s interval)
Memory	Yes
Additional heads	Clean-series
Brush dimensions	24 x 24 x 2,7 cm
Brush weight	119 g
Waterproof	IPX8
Input voltage	5 V – 0,5 A
Battery	2000 mAh, Li-ion
Charging time	4-5 hours

CASE SPECIFICATIONS

Battery life	17 days (2x daily 5 minutes = 34 sterilizations)
Charging	USB-C cable
UV LEDs	1x UV LED (UVC + UVA)
UV wavelength	265 – 280 nm + 395 nm
Sterilization time	5 minutes
Case dimensions	22 x 6,3 x 3,4 cm
Case weight	132 g
Input voltage	5 V – 350 mAh
Battery	600 mAh, Li-ion
Charging time	2-3 hours

MANUFACTURER'S LIABILITY

TrueLife products are warranted against manufacturing defects. Obligation of elem6 s.r.o. in the event of such defects, it is within the warranty to repair or replace the defective part or product with a comparable part or product at its own discretion. Remodeling, alteration or other unauthorized handling of the product means exceeding the scope of this warranty. This is even if the product is lost or damaged as a result of a defect, negligence or other error. The warranty does not cover damage resulting from use, accident or normal wear and tear. The company elem6 s.r.o. assumes no responsibility for any accident, injury, death, loss or other claim related to the product or resulting from the use of this product. The company elem6 s.r.o. shall in no event be liable for incidental or consequential damages related to or arising from the use of this product or any of its components. Returns or exchanges of spare parts and/or products may be subject to shipping, handling, exchange and/or restocking fees.

SAFETY WHEN WORKING WITH THE EQUIPMENT:

- Read the package safety information carefully.
- For proper handling of the product, read the current version of the manual, which is available for download at: www.truelife.eu
- Typographical errors and information changes reserved.

Manufacturer:
elem6 s.r.o.,
Bráskovská 15, 161 00 Prague 6, Czech Republic

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