EN – USER MANUAL

Before first use, the consumer must first read the instructions and safety precautions of the product.

1. Product description


Parts of the thermometer:

- Infrared sensor
- LCD display
- MODE / SET button
- UP button
- DOWN button
- BACKLIGHT
- SCAN button
- Battery cover

2. Display description

Display data:

- Surface mode (surface temperature measurement)
- Body temp mode (body temperature measurement)
- Temperature reading
- Memory number of the reading
- Beeper
- Temperature scale °C / °F
- Low battery warning
- Summary of stored readings

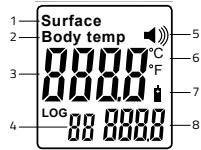3. Main technical specifications

Units of measurement	0,1 °C (0,1 °F)
Operating conditions – ambient temp.	10–40 °C (50–104 °F)
Storage conditions – ambient temp.	0–50 °C (32–122 °F)
Operating conditions – humidity	≤ 85%
Storage conditions – humidity	≤ 90%
Power source	2x 1.5V AAA batteries
Dimensions	98 x 36,5 x 132 mm
Weight	73 g
Measurement range – body temp mode	32,0–43 °C
Measurement range – surface temp mode	0–100 °C
Max. reading variation – body temp	+/- 0,3 °C
Max. reading variation – surface temp	+/- 1,0 °C
Distance from sensor to object	1–3 cm
Measurement speed	≤ 0,8 seconds
Automatic off	15 seconds

4. Accuracy of body temperature readings

32,0 – 35,9 °C (93,2 – 96,6 °F)	± 0,3 °C (32,0 – 34,9 °C)
36,0 – 39,0 °C (96,8 – 102,2 °F)	± 0,2 °C (35,0 – 42,0 °C)
39,0 – 42,5 °C (102,2 – 108,5 °F)	± 0,3 °C (42,1 – 43,0 °C)

5. Normal body temperature ranges by measurement method

Measurement method	Normal temperature
Anally	36,6 – 38,0 °C (97,8 – 100,4 °F)
Orally	35,5 – 37,5 °C (95,9 – 99,5 °F)
In the armpit	34,7 – 37,3 °C (94,4 – 99,1 °F)
In ear	35,8 – 38,0 °C (96,4 – 100,4 °F)

6. Normal body temperature ranges by age

Age	Normal temperature
0 – 2	36,4 – 38,0 °C (97,5 – 100,4 °F)
3 – 10	36,1 – 37,8 °C (97,0 – 100,0 °F)
11 – 65	35,9 – 37,6 °C (96,6 – 99,7 °F)
> 65	35,8 – 37,5 °C (96,4 – 99,5 °F)

7. Thermometer settings

- Setting measurement mode: Press MODE to toggle between surface temperature reading (it displays surface) and body temperature reading (it displays body temp). Press and hold MODE / SET for 2 seconds to access further settings. Press MODE / SET again to scroll through the individual settings.
- Pressing and holding MODE / SET for 2 seconds displays F-1: Press UP or DOWN to set the reading to °C or °F respectively.
- Pressing MODE / SET again displays F-2: Press UP and DOWN to adjust the temperature threshold for the beeper in increments of 0,1 °C.
- Pressing MODE / SET again displays F-3: Press UP or DOWN to enable or disable the elevated temperature beeper. 0' means OFF and 1' means ON.
- Pressing MODE / SET again clears the display and the settings are saved.


8. How to take accurate readings

Hold the thermometer 1–3 cm from the surface of the object or forehead. Press SCAN to take the reading. The measurement is displayed within 1 second.

Important notice:

- If the thermometer has not been used for a long time it will take 2 seconds longer to start up while the thermometer measures the current ambient temperature and calibrates.
- If there is sweat, cosmetics, hair or other obstructions on the forehead the thermometer may show measurement variations. In this case, take the reading from behind the ear, making sure there are no obstructions.
- If the ambient temperature is outside the operating range, it will display the message ERR: This means that any reading may have a greater variation than specified in the manual.
- If the body temperature or object temperature are outside the operating range it will display HI for higher and LO for lower than the operating temperature range.

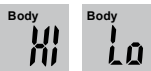9. Replacing the batteries

- The battery icon is on the display. If the battery icon flashes this means that the battery should be replaced. A low battery can cause inaccurate results. Open the battery cover, replace the battery and put the cover back on. Insert the battery, making sure the positive and negative are the right way round. If the battery is inserted the wrong way round it may cause damage to the thermometer.
- Remove the battery when the thermometer will not be in use for a long time.
- The thermometer can take 40,000 readings per battery.
- Dispose of the batteries at a collection point or specially designated place. Do not dispose with household waste, in a fire or in water.

10. Troubleshooting

- The reading is unexpected
 - The thermometer should be used under stable conditions. The reading may have been taken under frequently changing ambient conditions or when a very hot and cold object have been scanned in quick succession.
 - Note the correct distance of the scanner from the object. It must be 1–3 cm.
- Nothing happens when the buttons are pressed
 - Battery may be dead. Replace the battery.
 - The battery may be inserted the wrong way round. Check that the positive and negative are the correct way round.
 - The device may be damaged. Contact the service centre.
- If the problems persist after correctly following the instructions and procedures in this manual, contact your dealer or service centre. Never attempt to repair the device yourself. This product is covered by a 2-year warranty, valid from the date of purchase.

Warning: This product is intended to be used in households and similar applications.

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EC	REP	ALTHOME France, Z.A. 32390 MONTESTR UC/gers, France
----	-----	--

Shenzhen Calibur Industries Co., Ltd.
No. 11 Bldg., Shiguan Industrial Park, Gongming, Guangming New District, Shenzhen, China

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. You can determine a warranty service at the place of purchase. Manufacturer: Shenzhen Calibur Industries Co., Ltd., B46 Bldg., 2nd Ind.zone, Tanglang, Xili, Shenzhen, 518055 China. In case of technical problems or queries contact your dealer who will inform you about the following procedure. Products are covered by a two year warranty, unless otherwise stated.

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